
TESLA

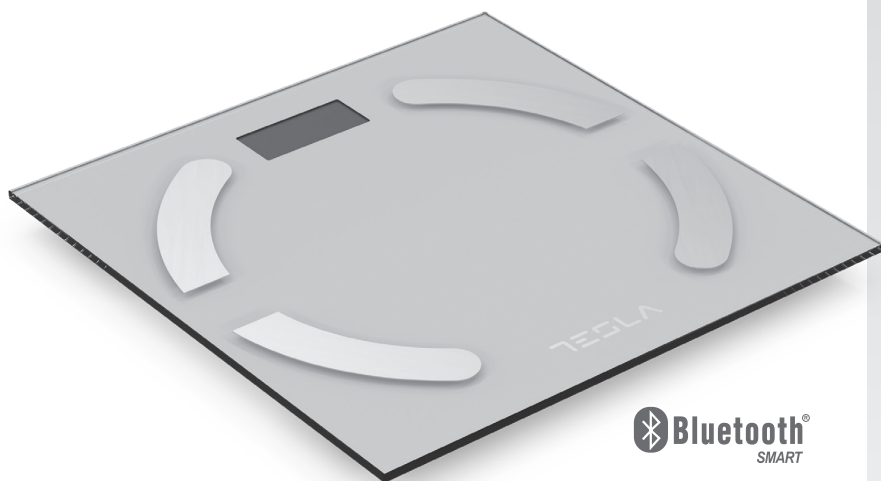BS301WX

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

Household Use only.
Read this instruction manual carefully before using.

Download now:
OKOK International

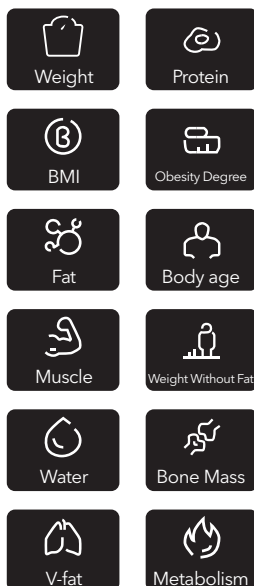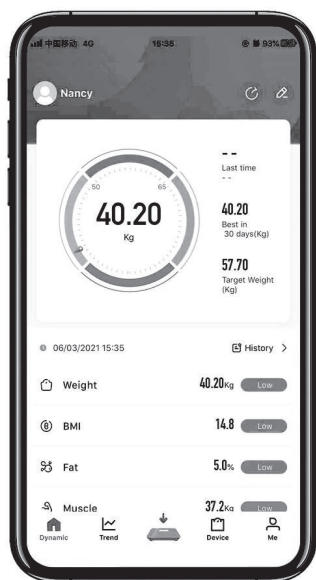HOW TO USE THE SCALE

The device needs to be connected with the app before first time use.

The weight scale will automatically connect to your phone and the application.

Be sure to activate Bluetooth on your phone or tablet before using the scale.

1. Connect your phone or tablet to the Internet.
2. Use a QR code scanner or search OKOK International from the App Store / Google Play to install APP.

Available on the
App Store

ANDROID APP ON
Google play

3. Pair the device, the scale will automatically connect to the phone and the application.
4. Register and complete user information (enter name, gender, height, birthday and desired weight) Click on the following to confirm your personal information.
5. Step on the scales and start measuring. Wait until the device shows your measurements.
6. Synchronized automatically, see the result.

- The device will display body weight, BMI, fat, muscle, water, visceral fat, bones mass and metabolism. The device will save your details automatically.
- If the Bluetooth is disconnected, data will be saved into the memory of the scale, sync to app while Bluetooth is connected.
- The user can type weight data into the app to keep a track

- The device will display current weight and generate reports and give health advice based on measurements. Including current weight interface, weight list interface, date report and notice.
- Weight trend will include weekly report, monthly report and quarterly report, indicating changes of weight through curves.
- Account management includes registration, login, changing password, forgot password, links to social media, and more.
- User management includes user information like nickname, birthday, height, gender, etc

Weighing Function

- Place the scale on a hard, flat surface (avoid carpets). An uneven floor contributes to inaccurate measurements.
- Before using the scale it is necessary to let it be calibrated and when the display turns off you can measure yourself.
- Get on the scale slowly. Put your feet on the same level without moving and wait for your weight on the screen.
- Distribute your weight equally and wait for the scale to determine your weight.
- The digits on the display will flash before your exact weight will be shown.
- After you get down the scale, it will be switch-off automatically after 30 seconds feet on the same level without moving and wait for your weight on the screen.

Cleaning the scale

- Do not immerse the balance in water: the electronic components can be damaged.
- Do not use a wire brush or any abrasive items.
- Do not use any acidic or abrasive detergents.
- Clean the device with a slightly moist cloth. Make sure that no humidity or waters seeps into the device!
- Use a commercially sold detergent or soap for cleaning.
- Dry the device well after cleaning with a dry and soft cloth.

Battery warning

- Do not mix old and new batteries.
- Do not mix alkaline, standard (carbon-zinc), or rechargeable (nickel-cadmium) batteries.
- Insert batteries using the correct polarity.
- Do not short-circuit the supply terminal.
- Remove batteries when not in use. Recycle or dispose of exhausted batteries according to local laws. Do not dispose of batteries in a fire - they may leak or explode.

Disposal

Do not Dispose of batteries with house-hold waste!

Batteries and packaging should be sorted for environmentally friendly recycling.

Only for EU countries:

according to the Directive 2006/66/EC, defective or used batteries, must be collected separately and disposed of in an environmentally correct manner.

CAUTION

- Handle the balance with care, as it is a precise measuring instrument. Do not drop it or jump on it.
- Do not load the appliance with more than 180 kg to avoid damaging the device.
- Do not leave your scale with water on it, it may damage the internal parts
- If the scales are carried or not going to be used, please store them only in a horizontal position.
- Use your scale with caution. This is a High precision strain gauge sensors system. Do not drop it and do not jump on it
- Do not get on it with wet feet, you could fall down
- The body fat analysis scale is not suitable for the following persons:
 - Children under 10 years and adults over 99 years
 - Persons with symptoms of fever, edema, or osteoporosis
 - Persons undergoing treatment for dialysis
 - Persons who use a pacemaker
 - Persons who take cardiovascular medication
 - Pregnant women
 - Sportspersons who perform over 5 hours of intensive sporting activities each day
 - Persons with a heartbeat lower than 60 beats/minute

Technical specifications

Product description: Bluetooth Body Fat

Scale Material: Tempered glass, plastic, Stainless steel

Product size: 300x300x28mm

Weight: 1.8 kg

Battery: 2x 1.5 V AAA Batteries

Capacity: Maximum weight 180 kg

Weight unit: Kg/lb

Memory function: Up to 8 persons

Analysis: Weight/ BMI/ Fat/ Muscle / Water/ V-fat/ Bone/ Metabolism / Protein....

Indicator: Overload Indicator, Low battery indication

App name: OkOk international

OKOK INTERNATIONAL APP

Support Android 4.3 and above for selected Android devices with Bluetooth 4.0 and above Install the app. Connect the unit with your mobile or any other device via Bluetooth. The analysis will be shown in the app called "OKOK" which you can easily download from the app store or google play.

The app can save data for up to 8 people, So you and your family and friends can use this scale

EU DECLARATION OF CONFORMITY

Manufacturer business name: TESLA/COMTRADE DISTRIBUTION D.O.O. BEOGRAD

Manufacturer address: BULEVAR ZORANA ĐINĐIĆA 125i, 11000 BEOGRAD

WE, MANUFACTURER OF BELOW DESCRIBED PRODUCTS, HEREBY DECLARE UNDER OUR SOLE RESPONSIBILITY, THAT THE PRODUCTS TO WHICH THIS DECLARATION RELATES:

product description:

Brand	Model	Type of products
TESLA	BS102B, BS103W, BS301WX	Body Scales

ARE IN CONFORMITY WITH THE FOLLOWING STANDARDS AND LEGISLATIVE ACTS:

Harmonized standard	Legislative act
EN 55014-1:2017+A11:2020 EN 55014-2:2015	<i>Electromagnetic Compatibility (EMC)</i> 2014/30/EU

Information on authorized representative in the EU who has a written mandate from the above stated manufacturer to act on its behalf in carrying out certain tasks required in the applicable Union harmonization legislation including tasks specified under Article 4 of Regulation (EU) 2019/1020:

Business name of representative in EU:

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E-mail address of representative in EU:

office@tesla.info

Place and date of issue: *Belgrade, 23.08.2023*

Signed for and on behalf of manufacturer:

(Signature and stamp)

Name of authorized person:

Nebojša Lozo

Function of authorized person:

General manager

TESLA

tesla.info