TESLA

IONIC POWERED
TECHNOLOGY

PREMIUM Home Toothbrush
TDP121PW, TDP121NB

User Manual

ENG

ALB

BIH/CG

BG

EST

GR

HR

HU

LT

LV

MK

RO

SLO

SRB

ELECTRIC TOOTHBRUSH

Read this instruction manual carefully before using.

Thank you for choosing Tesla.

We are committed to providing highly cost-effective products to help clean teeth and protect oral health. This brush can effectively reduce the accumulation of plaque, coffee, tea, smoking stains, and more. It can also protect gums and play an important role in improving and preventing gingivitis.


IMPORTANT SAFEGUARDS

- Using the Tesla Brush is a new kind of brushing experience.
- Do not open, repair or alter the product without authorization.
- Do not strike or drop the product. If product appears defective, please stop using it.
- Please check whether local voltage matches the designated voltage on the power supply before connecting.
- Use this product only with the designated Tesla toothbrush head and charger.
- If you've had past oral or gum surgery, or are experiencing any sort of oral discomfort, please consult with your dentist before use.
- Though use of the Tesla Home toothbrush is generally safe for all, if you wear a cardiac pacemaker or similar device, please consult your physician before using this brush.
- This product was designed to clean only teeth, gums and tongue.
- To ensure proper usage, persons (including children) inexperienced with the use of this type of product should not use it without supervision.
- Do not use the charging base outdoors or under high temperature. Do not immerse the base or body in water or any liquid. Please remove power supply when not in use and before cleaning it. Please do not clean with hair care products (such as gels or oils), hand soap, or acidic cleaning fluid, as these will damage the surface of the product.
- Children should be supervised to ensure proper usage, and that they do not play with or mishandle the toothbrush.
- Do not store under high temperature conditions or direct sunlight.

PRODUCT INFORMATION

Product name	Tesla Home Toothbrush	Model	TDP121PW, TDP121NB
Handle power	1.2 W	Handle working voltage	DC 3.7 V
Charging power	1.5 W	Charger working voltage	DC 5 V
Charging time	Approximately 4 hours	Battery capacity	260 mAh
Battery usage on full charge	100 minutes	Weight	Approximately 85g (Including brush)
Body size	Approximately 2,54 cm x 2,54 cm x 22,86 cm		
Frequency	Approximately 22.000 strokes/min. (Clean mode); Approximately 26.000 strokes/min. (White+ mode)		

COMPONENT NAME

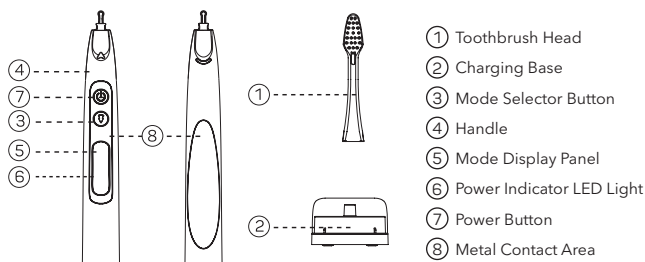
USING THE TESLA HOME TOOTHBRUSH

1. Install the replaceable brush head firmly over the steel shaft of the unit's handle.
2. With your moisturized hand, hold the handle, while making contact with the metal plate. Doing this is what generates the negative ion flow that makes IONPA sonic power work more effectively.
3. Wet the bristles and apply a small amount of toothpaste, if desired. (If toothpaste is used, less than half your regular amount of toothpaste is generally sufficient.)
4. Gently place the bristles of the brush on your teeth, angled toward the gum line.
5. Select the appropriate mode by pressing the Mode Selector button. The display panel will show which mode is selected. Choose from four modes, as follows:
 - a. Clean:** approx. 22,000 strokes per minute. This is an everyday mode that addresses all matters of oral hygiene.

- b. White+:** approx. 26,000 strokes per minute, for when you wish to give special attention to eliminating surface stains.
 - c. Gum care:** In addition to polishing teeth, this mode provides gentle, massaging stimulation to the gums.
 - d. Hi-Clean:** This mode is slightly more vigorous, to get at plaque, deep in the recesses between teeth.
6. Gently move the brush head slowly across the teeth in a slight back and forth motion to achieve the best clean. Alternate brush strokes between front, back, and chewing surfaces of teeth, devoting at least 10 seconds to each surface. Be sure to allow the bristles to reach between teeth and massage gums. Do not place too much pressure on the bristles against the teeth, allow the strength of the oscillation to do the work.

Replacement brush head indicator


If you choose not to use replacement indicator:

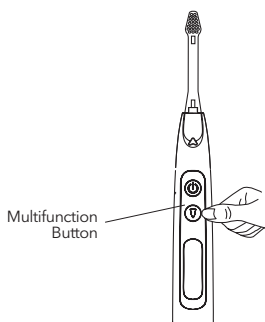The replacement indicator is active when you open the package. If you choose not to use the indicator, first make sure your device is off and press the multifunction button for 2 seconds to turn off the indicator.

How to reset the replacement indicator

(Make the brush head icon disappear)

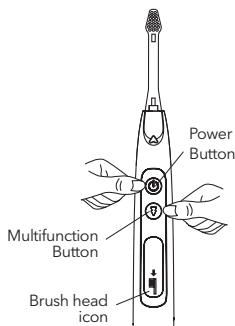When the indicator is active, the brush head icon will appear after 180 minutes of use. The icon will blink for 30 seconds even after the device is turned off.

To reset, turn the device off and press both the power button and the multifunction button together for 2 seconds. The indicator will not reset by simply replacing a brush head.

FEATURES


1. The Tesla motor delivers up to approx. 26,000 strokes per minute, with an amplitude (coverage sweep) of 6 mm, to clean teeth effectively.
 2. The brush's body is ergonomically designed with a soft handle finish which gives the toothbrush a comfortable feel and a reliable grip.
 3. The Smart Timer pauses every 30 seconds, three times in every two-minute session, to remind the user of the recommended two-minute optimal brushing time.
 4. The body of the brush is immersible to one meter for 30 minutes without leakage, more than sufficient to meet daily cleaning requirements.
 5. The unit contains a safe and durable lithium battery. The fully charged battery can run for 25 days (two-minute usage twice per day) and alerts the user when the charge is low.
 6. It offers 4 modes, Clean, White+, Gum care, and Hi-Clean, to meet different needs.
- The unit's functions were factory pre-set and checked before distribution.

INSTRUCTIONS FOR REPLACING THE BRUSH HEAD

DETACH

Twist the head 90 degrees as shown in the image on the left and then pull apart from the handle.

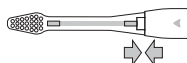INSTALL

Notch on the neck should connect with inner notches of toothbrush head. Position new brush head and slide down the rod until it clicks into position.

CHARGING THE UNIT

1. Thanks to inductive charging, simply placing the toothbrush on the AC-connected charging base will maintain the battery charge.
2. LED light indicates that the unit is charging. This indicator light will go off when the unit is fully charged, which should require no more than four hours.
3. The Mode Selector button is disabled during charging.
4. If you place the active device on the charging base, it will stop working and enter charging mode.

CLEANING THE UNIT

1. Handle

Remove the brush head from the handle. Use a soft, wet cloth to wipe the entire handle and shaft surface and allow to air dry.

2. Charging base

Use a soft, wet cloth to wipe the charging base and air dry.

- Do not immerse in water.

3. Brush heads

- Fully rinse brush head and bristles after each use.
- Avoid squeezing or pull the bristles, to avoid bending them out of shape.
- We recommend changing the toothbrush head every 45 days. Overused brushes can accumulate bacteria, and the cleaning effectiveness will diminish greatly over time.

Dispose responsibly in accordance with local regulations

The battery must be removed from the product before disposal. Removed batteries must not be disposed of with household waste. Each consumer is legally obliged to dispose of batteries in accordance with national regulations. Keep removed batteries out of reach of children; do not throw into fire, short-circuit, or disassemble.

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