TESLA

IONIC POWERED MANUAL TOOTHBRUSH

Tesla Kids Brush

TMT8

TMT8P

TMT8Y

TMT9

TMT9P

TMT9Y

User Manual

ENG

ALB

BIH/CG

BG

EST

GR

HR

HU

LT

LV

MK

RO

SLO

SRB

IMPORTANT SAFEGUARDS


- Parents: Teach younger children the proper use of this toothbrush.
- This product is designed for home use to clean teeth, gums, and tongue. Do not use it for any other purpose.
- Do not open, repair, or modify the product without authorization.
- Do not soak the product in water for extended periods.
- Do not leave the product in high-temperature environments or expose it to direct sunlight.
- Children and persons with reduced physical, sensory, or mental capabilities, as well as those lacking sufficient experience and/or knowledge, may use the product only under supervision or after receiving instruction on safe use and understanding of potential risks.
- Do not allow children to play with the device.
- Discontinue use of the brush head if the bristles are broken or bent. Replace the brush head as soon as signs of wear appear.
- Do not wash any parts of the device in a dishwasher.
- If battery leakage occurs, avoid contact with the liquid on skin, eyes, or mucous membranes. In case of contact, rinse the affected area with plenty of water immediately and seek medical advice.
- Consult your dentist before using the device if you have had oral or gum surgery in the past two months.
- If you have a pacemaker or any other implant, consult your doctor or the implant manufacturer before using the device.

Dispose responsibly in accordance with local regulations

The battery must be removed from the product before disposal. Removed batteries must not be disposed of with household waste. Each consumer is legally obligated to dispose of batteries in accordance with national regulations. Keep removed batteries out of reach of children; do not throw into fire, short-circuit, or disassemble.

TOOTHBRUSH PARTS

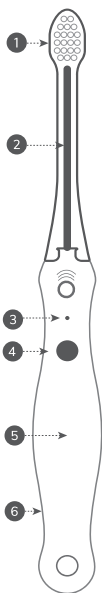- 1 Replacement head
- 2 Metal rod for transfer of ions
- 3 Green light battery indicator
- 4 Battery status check button
- 5 3-volt battery in a waterproof compartment covered with metal tape
- 6 Black TPE stripe

Clinically proven

Tesla Kids Brush removes dental plaque up to 48% more efficiently than other manual toothbrushes.

48%

Replacement heads

Change your replacement head every 2 or 3 months for the best result.

3 Volt battery

Lithium battery is integrated in the body of the toothbrush, and it lasts more than one year.

Japanese Technology

Synonym for quality and innovation. Tesla Kids Brush, officially certified by Tokyo Institute of Technology.

Ionization Technology

Negatively charged ions successfully loosen the bond between teeth and dental plaque.


No toothpaste required

First toothbrush that can be used efficiently without toothpaste.

HOW IT WORKS?

Ionization

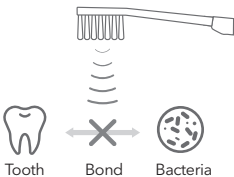How does Tesla Kids Brush remove dental plaque?


Dental plaque is attached to teeth with positive ions from saliva. Negative ions from the Tesla Kids Brush effectively loosen these bonds between teeth and plaque and allow easy removal of plaque by brushing.

How to use it?

Tesla Kids Brush is very simple to use.

While you hold your toothbrush, keep contact with the black TPE stripe with your moist finger or palm.

And brush your teeth properly.

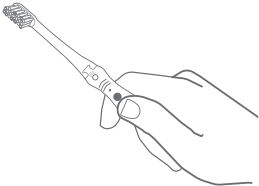How to brush your teeth properly?

Always repeat the same action.

Brush outside first, then all of the inside. Do the lower parts first, and then the upper.

On the outside, place the toothbrush bristles against the teeth at a slight angle toward the gum line and brush in gentle short circular movements. Take 5-10 moves for every tooth.

For the inside parts, use the same technique, moving from one side molar to the other.

As for chewing sides, start at the one side molar and work your way forward slowly in short back-and-forth movements to the other side molar.


Whole process should take between 3 and 6 minutes.

Clinically proven, Tesla Kids Brush removes dental plaque up to 48% more efficiently than other manual toothbrushes.

There is no need to apply any pressure. Brush gently.

Tesla Kids Brush is clinically proven to be effective against dental plaque even without toothpaste. Please advise your dentist regarding usage of toothpaste.


How to check the battery condition?

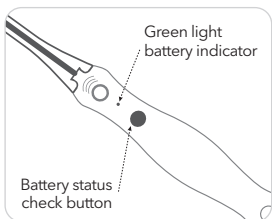When you want to check if the battery is sufficiently charged, just press the battery status check button.


Green light indicates that your toothbrush is ready for use. Lithium coin cell battery (3V) is integrated in the body of the toothbrush and it lasts more than one year. When battery is depleted dispose responsibly according to local regulation.

How to remove the replacement head?

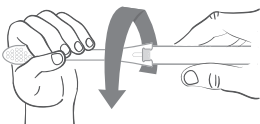Hold the handle firmly. Rotate the replacement head 90 degrees and pull it.

How to attach the replacement head?

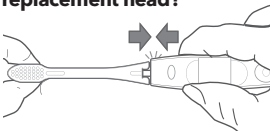Slide the replacement head onto the metal rod until it snaps into position.

You can find different types of Tesla Brush replacement heads. Check out www.tesla.info for further details.

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