TESLA

IONIC POWERED MANUAL TOOTHBRUSH

Tesla Original Brush

TMT1
TMT1G
TMT1P
TMT10

TMT2
TMT2P
TMT2W
TMT3

TMT3P
TMT3W
TMT4
TMT4P

TMT5
TMT5P
TMT6
TMT7

User Manual

ENG

ALB

BIH/CG

BG

EST

GR

HR

HU

LT

LV

MK

RO

SLO

SRB

IMPORTANT SAFETY MEASURES


- Before first use, carefully read this manual and follow all safety instructions.
- This device is intended exclusively for cleaning teeth, gums, and tongue.
- Consult your dentist before using the device if you have undergone oral or gum surgery in the past two months.
- If you have a pacemaker or another type of implant, consult your physician or the implant manufacturer before using the device.
- Use only original accessories and consumable parts from Tesla.
- Children and persons with reduced physical, sensory, or mental capabilities, as well as those lacking sufficient experience and/or knowledge, may use the product only under supervision or after receiving instruction on safe use and understanding the potential risks.
- Children are not allowed to clean or maintain the device on their own and may do so only under adult supervision.
- Do not allow children to play with the device.
- Stop using the brush head if the bristles are broken or bent. Replace the brush head at the first signs of wear.
- Do not wash any parts of the device in a dishwasher.
- Never attempt to recharge batteries that are not designed to be rechargeable.
- Keep the device and batteries away from open flames, direct sunlight, and high temperatures.
- If the battery leaks, avoid contact with the fluid on skin, eyes, or mucous membranes. In case of contact, immediately rinse the affected area with plenty of water and seek medical advice.
- Do not place the device or batteries in microwave ovens or on induction cooktops.
- Do not open, repair, or modify the device without proper authorization.

Dispose responsibly in accordance with local regulations

The battery must be removed from the product before disposal. Removed batteries must not be disposed of with household waste. Each consumer is legally obliged to dispose of batteries in accordance with national regulations. Keep removed batteries out of reach of children; do not throw into fire, short-circuit, or disassemble.

TOOTHBRUSH PARTS

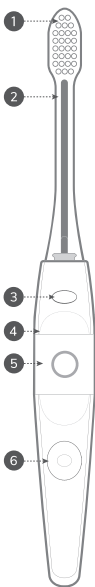- 1 Replacement head
- 2 Metal rod for transfer of ions
- 3 Red light battery indicator
- 4 3-volt battery in a waterproof compartment covered with metal tape
- 5 Metal tape
- 6 Battery status check button

Clinically proven

Tesla Brush removes dental plaque up to 48% more efficiently than other manual toothbrushes.

48%

Replacement heads

Change your replacement head every 2 or 3 months for the best result.

3 Volt battery

Lithium battery is integrated in the body of the toothbrush, and it lasts more than one year.

Japanese Technology

Synonym for quality and innovation. Tesla Brush, officially certified by Tokyo Institute of Technology.

Ionization Technology

Negatively charged ions successfully loosen the bond between teeth and dental plaque.


No toothpaste required

First toothbrush that can be used efficiently without toothpaste.

HOW IT WORKS?

Ionization

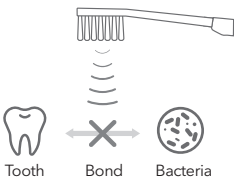
How does Tesla Brush remove dental plaque?

Dental plaque is attached to teeth with positive ions from saliva. Negative ions from the Tesla Brush effectively loosen these bonds between teeth and plaque and allow easy removal of plaque by brushing.

How to use it?

Tesla Brush is very simple to use. While you hold your toothbrush, keep contact with the metal strip with your moist finger or palm. And brush your teeth properly.

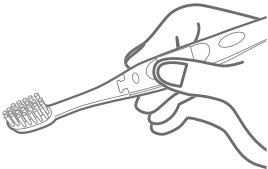How to brush your teeth properly?

Always repeat the same action. Brush outside first, then all of the inside. Do the lower parts first, and then the upper.

On the outside, place the toothbrush bristles against the teeth at a slight angle toward the gum line and brush in gentle short circular movements. Take 5-10 moves for every tooth.

For the inside parts, use the same technique, moving from one side molar to the other.

As for chewing sides, start at the one side molar and work your way forward slowly in short back-and-forth movements to the other side molar.


Brush for about 2 minutes for optimal results.

Clinically proven, Tesla Brush removes dental plaque up to 48% more efficiently than other manual toothbrushes.

There is no need to apply any pressure. Brush gently.

Tesla Brush is clinically proven to be effective against dental plaque even without toothpaste. Please advise your dentist regarding usage of toothpaste.


How to check the battery condition?

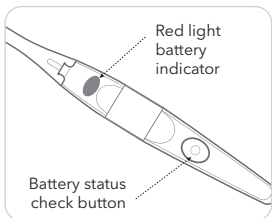When you want to check if the battery is sufficiently charged, just press the battery status check button.


Red light indicates that your toothbrush is ready for use. Lithium coin cell battery (3V) is integrated in the body of the toothbrush and it lasts more than one year. When battery is depleted dispose responsibly according to local regulation.

How to remove the replacement head?

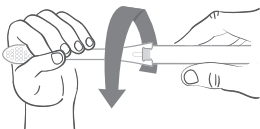Hold the handle firmly. Rotate the replacement head 90 degrees and pull it.

How to attach the replacement head?

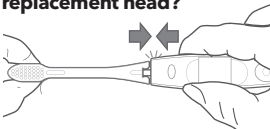Slide the replacement head onto the metal rod until it snaps into position.

You can find different types of Tesla Brush replacement heads. Check out www.tesla.info for further details.

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