TESLA

ELECTRIC TOOTHBRUSH

Tesla Travel

TDM011PW TDM011PG
TDM011RD TDM011BK

User Manual

ENG

ALB

BIH/CG

BG

EST

GR

HR

HU

LT


LV

MK

RO

SLO


SRB

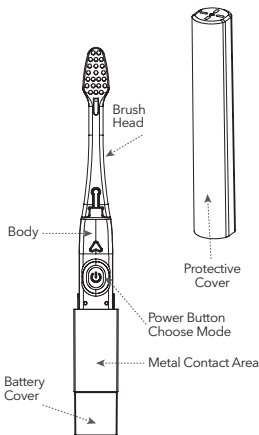
INSTRUCTIONS

Replacing the battery

Turn battery cover as shown in 1 and pull-down as shown in 2. Insert one new AAA battery, negative side first. Replace battery cover as shown (3 & 4).

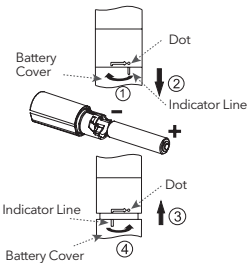
Changing the replaceable brush head

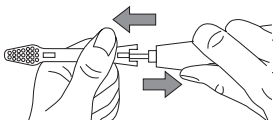DETACH

Grasp the brush head firmly and pull it straight off.

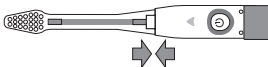INSTALL

Snap a new brush head into place over the silver bar.

Brushing Instructions

Power
Button

Press the power button to start the toothbrush; press again within 3 seconds to enter the massage function; pressing a third time, or after 3 seconds, will turn the toothbrush off.

With your hand slightly wet, hold the handle, while making contact with the metal plate. Doing this is what generates the negative ion flow that makes Tesla sonic power work more effectively.

Wet the bristles and apply a small amount of toothpaste, if desired. (If toothpaste is used, less than half your regular amount of toothpaste is generally sufficient).

Divide oral area into upper, lower, right and left sections.

Place the toothbrush bristles against the teeth at a slight angle toward the gum line.

Press the power button to turn on the toothbrush, gently move the brush head across the teeth with a small up and down movements. When you depress the Power Button the device automatically starts in clean mode. This is an everyday mode that is like using any manual toothbrush, but with the added advantage of sonic power. Press the Power Button once again, and the device enters Gum care mode. Use the brush to gently massage and stimulate the gums while cleaning teeth.

Timer function: The brush will pause every 30 seconds, three times every two-minute session, to remind the user of the recommended brushing time.

Brush each section for 30 seconds, with 10 seconds each on front, back and chewing surface.

To avoid splashing, turn toothbrush off before removing from the mouth by pressing the on/off button.

Clean the head and body of the brush after every use.

CAUTION!

- Parents: Teach younger children the proper use of this toothbrush.
- This product is designed for home use to clean teeth, gums, and tongue. Do not use it for any other purpose.
- Do not open, repair, or modify the product without authorization.
- Do not soak the product in water for extended periods.
- Do not leave the product in high-temperature environments or expose it to direct sunlight.
- Children and persons with reduced physical, sensory, or mental capabilities, as well as those lacking sufficient experience and/or knowledge, may use the product only under supervision or after receiving instruction on safe use and understanding of potential risks.

- Do not allow children to play with the device.
- Discontinue use of the brush head if the bristles are broken or bent. Replace the brush head as soon as signs of wear appear.
- Do not wash any parts of the device in a dishwasher.
- If battery leakage occurs, avoid contact with the liquid on skin, eyes, or mucous membranes. In case of contact, rinse the affected area with plenty of water immediately and seek medical advice.
- Consult your dentist before using the device if you have had oral or gum surgery in the past two months.
- If you have a pacemaker or any other implant, consult your doctor or the implant manufacturer before using the device.

PRODUCT SPECIFICATIONS

Power supply: 1 AAA alkaline battery

Battery life: 120 days (twice a day, two minutes per use)

Vibration:

approx. 22,000 stokes/min.

Size:

16.9cm (H) x 1,9cm (W) x 2,0 cm (D)

Weight:


approx. 50g (including the battery and brush head)

Waterproof grade: IPX7

DISPOSE RESPONSIBLY IN ACCORDANCE WITH LOCAL REGULATIONS

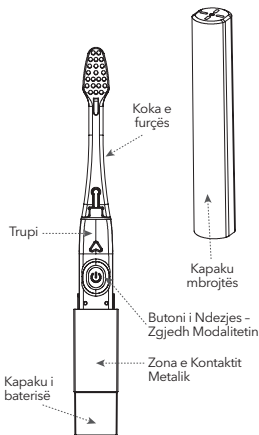The battery must be removed from the product before disposal. Removed batteries must not be disposed of with household waste. Each consumer is legally obligated to dispose of batteries in accordance with national regulations. Keep removed batteries out of reach of children; do not throw into fire, short-circuit, or disassemble.

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